



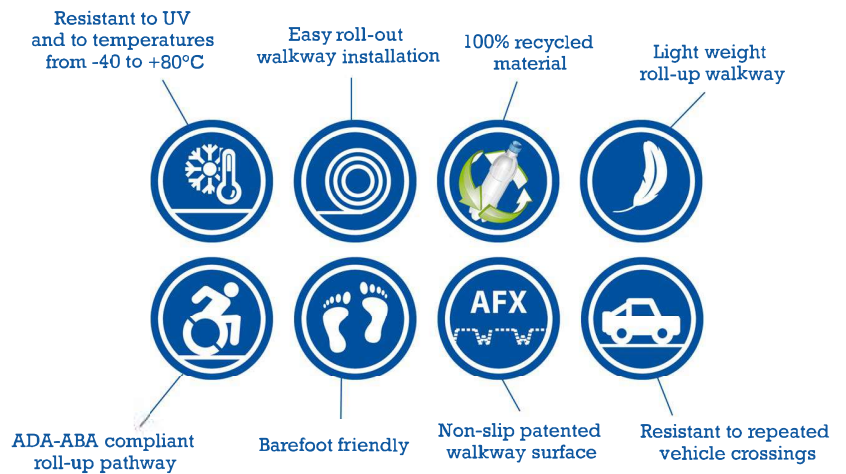

The Mobi-Mat is a portable and removable rollout access pathway for pedestrians and people who use wheelchairs, strollers, bicycles, mobility scooters and walking frames.

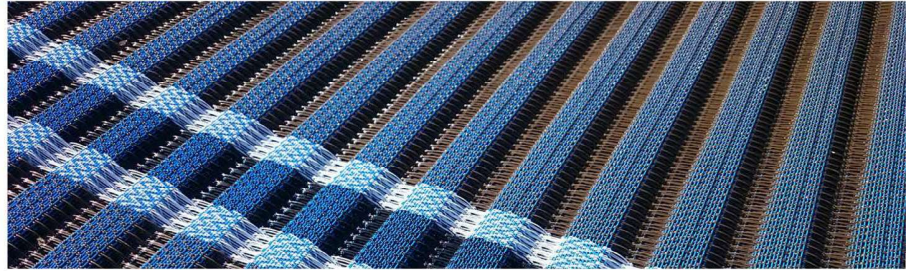

Mobi Mat provides a firm, safe, stable & smooth continuous surface without gaps.

Environmentally friendly - Made of 100% recycled polyester, Mobi-Mat contributes to the environment by utilising recycled plastics.

User friendly cool and smooth surface - The Blue Mobi-Mat remains cool when exposed to very high temperatures up to 80°C (176°F), for the benefit of all beach-goers.

Lightweight and Easy Installation - Mobi-Mat is 20.5kg for a 10m x 1.53m roll, making it easy to roll out with one or two people depending on the length of mat. Installation and retrieval can be performed by one person when using the Mobi-Roller.

Colour - Easily provides dedicated traffic lanes and is great for those with visual impairments.

Low Maintenance - The permeable structure allows sand to filter through without sinking into the sand. It is easily maintained by sweeping any excess sand build up with a broom or leaf blower.



Each kit is equipped with end connectors and eyelets to insert the staples provided in the kit.



Colour: Blue

Width: 1.53m

Length: 10m

Roll Weight: 20.5kg

Mobi-chair

BEACH ACCESS CHAIR

The Mobi-Chair is a high quality, floating beach wheelchair that provides beach and water access to people with limited mobility. The armrests act as flotation devices as you enter the water, providing the user with comfort and stability. It features quick-release components for easy transport and assembly.



Features:

Comfortable and Safe - Stable and well balanced with rubber pneumatic tyres absorbing shock and discomfort. Seatbelt for added safety.

Seat to Floor Height - Has a seat height of 19" which assists with transfers.

Durable - Manufactured with stainless steel and aluminium. Backrest fabric is non-allergenic, water and UV resistant. Waterproof back pocket.

Provides floatability - With the floatable armrests and tyres, Mobi-Chair easily transports the user from the beach to the water. Easy to push across hard and soft sand.

Assembly - Takes minutes to assemble and there are no tools required. Folds flat for easy storage.

Backrest recline - Recline adjustments set to 3 different positions, making it ideal for users with limited trunk control.



Dimensions:

Seat - 18" wide x 17" deep x 19" high

Chair - 60" long x 42" wide x 50" deep

Overall Weight -30KG

Front Castor -12" x 6.5" (Rubber Pneumatic)

Rear Wheel - 16" x 7.5" (Rubber Pneumatic)

Weight Limit - 136KG

