




Treat-Eezi Pressure Injury Overlays

Treat-Eezi works by using a unique 4-layer design to prevent and heal pressure ulcers, that also addresses skin shearing. Unlike standard overlays, it uses 4 layers of a three-way weave spacer fabric. This fabric evenly distributes pressure across the surface, while the layers move against each other to reduce skin shearing, assisting in the healing and lowering the risk of pressure ulcers.

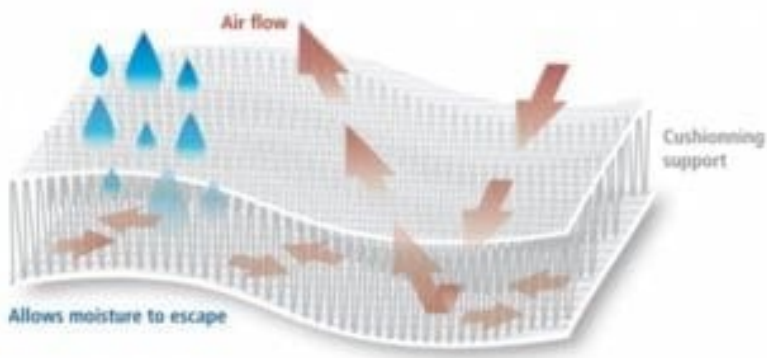

- Helps heal Stage 2, 3 & 4 pressure ulcers
- Pressure mapped at just 34mmHg
- Reduces shearing and minimises friction between skin & pad
- Helps regulate body temperature & wicks away moisture
- No electricity, no pumps, compact & easily transportable
- Maintenance free, can be washed to 92⁰ Celsius, can be bleached
- Allows sharing your bed with your partner
- Cost effective, easy to use - requires no training
- 159kg weight limit




C1 SOUTH

Treat-Eezi Mattress Overlay


Fits single beds and hospital beds.

195cm x 90cm

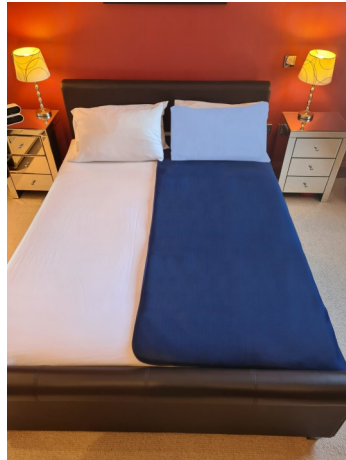

Treat-Eezi Small Overlay

Targeted skin protection without the need for a full overlay. 100cm x 80cm

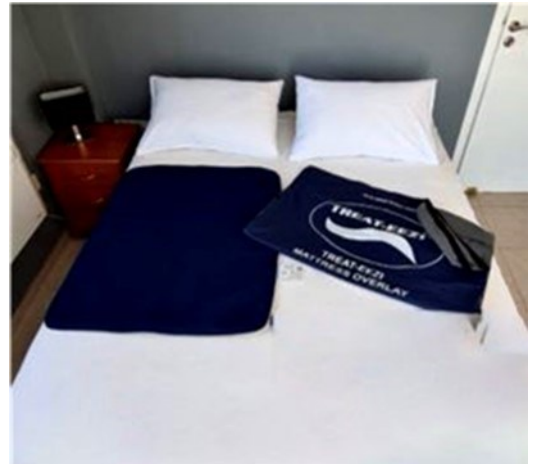

Treat-Eezi Heel Rest

Use on bed or on wheelchair footrest.

70cm x 30cm



Treat-Eezi Pillow Case

Fits up to standard European pillows size



Treat-Eezi Recliner Overlay

195cm x 60cm



Treat-Eezi Seat Pad

Wheelchairs, dining/lounge chairs or car.

40cm x 50cm



"We have used the Treat-Eezi Pressure Care Overlay in our nursing home for at least 3 months. We have found it to be highly effective in providing pressure relief with the added benefit of no noisy pump and no need for a contingency plan should we have a power cut." Manager - Park View Nursing Home - UK

Re: Treat-Eezi - "For the first time in over 10 years, he has started to have healing in his grade 4 pressure wound. Thank you!" Sue - OT - NZ

C1 SOUTH

0800 212 333 email: info@c1south.co.nz www.c1south.co.nz